

Cement-Retained Implant Restorations Made Easy with the STAGE-1 COC Quick-Pack

The STAGE-1 COC Quick-Pack delivers the required components for fabricating a cement-retained restoration in one convenient package. **The Quick-Pack includes:**



COC Abutment
Qty. 1



COC Quick-Cap
Impression Cap, Qty. 2



COC Abutment
Analog, Qty. 1




COC Plastic Sleeve,
Locking, Qty. 1


System Benefits:

- Three abutment heights virtually eliminate the need for abutment prepping
- Unique, one-piece impression cap retracts tissue and takes a precise impression in one easy step
- Color-coded components are simple to identify
- Resealable box ensures easy component transfer among the restorative team members
- Ideal for referral partnerships
- Components also available individually


Simple, accurate, reliable restorations:

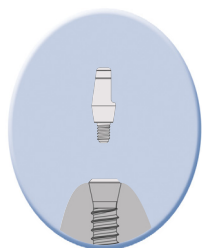


Seat Abutment

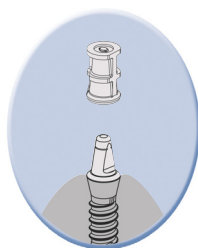


Place Quick-Cap

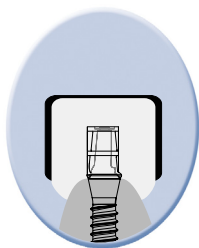

Take Impression

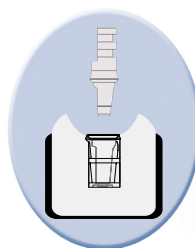

Insert Abutment Analog
into Quick-Cap

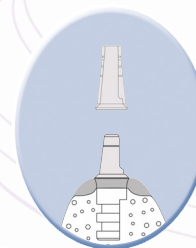

Place Plastic Sleeve on
Abutment Analog

Quick-Pack Ordering Information:

Regular Diameter (RDS) Prosthetic Connection

4.0mm Height Quick-Pack (RDS)	S2442-40
5.5mm Height Quick-Pack (RDS)	S2442-55
7.0mm Height Quick-Pack (RDS)	S2442-70

Wide Diameter (WDS) Prosthetic Connection

4.0mm Height Quick-Pack (WDS)	S2507-40
5.5mm Height Quick-Pack (WDS)	S2507-55
7.0mm Height Quick-Pack (WDS)	S2507-70



Manufacturer:

Lifecore Biomedical, Inc.
3515 Lyman Blvd.
Chaska, MN 55318 USA
Call: 800-752-2663
Fax: 800-651-8521

Outside the USA

Call: +1-952-368-4300
Fax: +1-952-368-4324

Authorized Representative:

Lifecore Biomedical, SpA
via A. Fleming, 19
37135 Verona, Italy
Call: +39-045-8230294
Fax: +39-045-8250296

Subsidiary:

Lifecore Biomedical, GmbH
Alfter, Germany
Call: +49-2222-9294-0
Fax: +49-2222-977356

Subsidiary:

Lifecore Biomedical, AB
Särö, Sweden
Call: +46-31-93-68-23
Fax: +46-31-93-68-45